

*Nasco* **Life/form**®

**Advanced “Airway Larry”  
Airway Management  
Trainer Head  
LF03684U  
Instruction Manual**



ISO 9001 and  
ISO 13485

*Nasco* Plastics  
ISO 9001 and ISO 13485

Products are manufactured using a quality management system  
certified as being in conformity with ISO 9001 and ISO 13485.

**Life/form**®

**5**  
Year  
Warranty

**Life/form**® Products by Nasco



## About the Simulator

The **Life/form**® Advanced Airway Management Trainer Head has the added features of laryngospasm and tongue swelling, this trainer will provide students with an easy-to-use model that will build their confidence and expertise while still providing a realistic intubation experience.

Landmarks include nostrils, teeth, tongue, oral and nasal pharynx, larynx, cricoid ring, epiglottis, arytenoid, false and true vocal cords, trachea, esophagus, “Airway Larry” lung set, and stomach.

Nasco’s Advanced “Airway Larry” allows you to practice oral, digital, and nasal intubation. E.T., E.O.A., PTL®, LMA, and Combitube® insertion can all be practiced as well (please see section titled “Using the Combitube®”). Suction techniques and proper cuff inflation can also be performed and evaluated.

In addition, this model also allows students to apply cricoid pressure during the intubation procedure (Sellick’s Maneuver). By locating the larynx and depressing the cricoid ring, the student may realistically lower the vocal cords farther into the field of view. This also effectively closes off the esophagus of the simulator.

The Airway Management Trainer head can be placed on **CPARLENE**® brand manikins creating a valuable mega-code component. With proper care, your **Life/form**® Airway Management Trainer will provide years of valuable service. Please review the instructions carefully for setup, cleanup, and maintenance procedures.

## List of Components

- Advanced “Airway Larry” Airway Management Trainer
- 8-oz. Pump Spray Lubricant
- Left and right simulated lungs
- Simulated stomach
- 20 cc syringe (2)
- 3 Red Caps

## Setup

The Advanced Airway Management Trainer is ready to use upon delivery. Simply remove it from the carton, and unwrap the packaging material. The lung and stomach bags are already in place and ready to use.

**Figure 1**



**Figure 2**




## Lubrication


Before using the Advanced Airway Management Trainer, lubricate both the simulator and supplies being used with the pump spray lubricant provided. **(See figures 1 & 2.)**



**Figure 3**

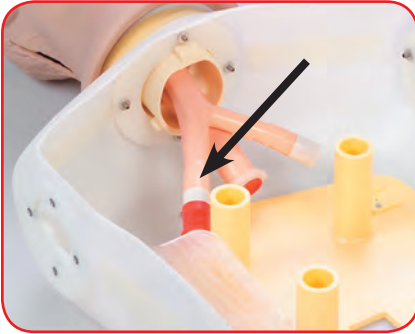

**Figure 4**

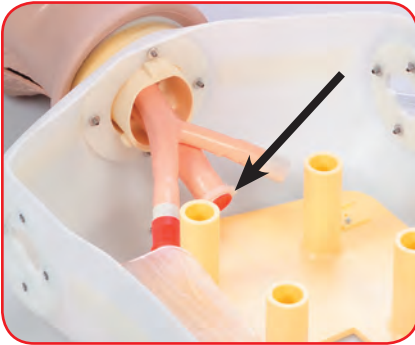

**Figure 5**

### Using the Combitube®

Thoroughly read and follow the instructions that come with the Combitube®. The trainer will accept either a full-size or a small adult tube. As with a live patient, it may be necessary to back the tube out slightly if ventilation cannot be established.

**Note:** Depending on tube placement, the large cuff may not accept the recommended amount of air. In this case, simply inflate the cuff to its maximum volume (when the plunger stops), detach the syringe from the blue pilot balloon, and proceed.

### Attaching the Head to the Manikin

As explained earlier, Nasco's Advanced Airway Management Trainer has the ability to be attached to **CPARLENE**® and modular type manikins. To attach the Airway Management Trainer Head to the above manikins, you must first remove the lung and stomach bags if they are attached to the trainer.

Next, turn the head face down and begin inserting it at the shoulders. Tip the head upward slightly until the two keyways on the back of the neck slide under the ring and into the top of the torso. **(See figure 3.)**

Now match the single keyway on the neck front up to the slot on the backside of the torso and push the head down until it snaps into place. Rotate the head back around so that it is facing up and in the correct position.

### Connecting the Airways

Once the head is in place, connect the **CPARLENE**® lung bag to the right bronchus by sliding the cap plug over the white adapter in the bronchus. **(See figure 4.)**

Then use the provided red cap plugs to plug the esophagus and the left bronchus. **(See figure 5.)** This will ensure the system is sealed and operational.

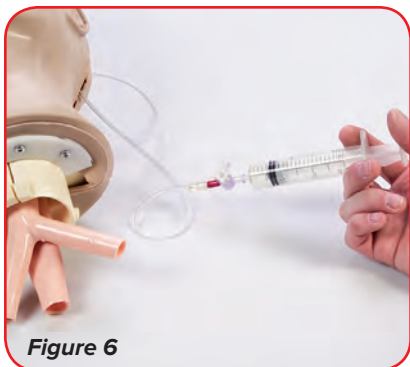


Figure 6

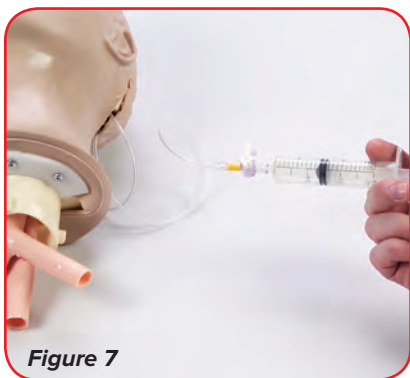


Figure 7

## Laryngospasm

Pull back the plunger of a 20 cc syringe and attach it to the red stopcock. **(See figure 6.)** Depress the plunger to force 1-5 cc of air into the system. Close the tubing port at the stopcock to maintain the pressure within the larynx. The syringe may be removed at this point.

## Tongue Swelling

Pull back the plunger of a 20 cc syringe and attach it to the yellow stopcock. **(See figure 7.)** Depress the plunger to force approximately 10 cc of air into the system. Close the tubing port at the stopcock to maintain the pressure. The syringe may be removed at this point.



Figure 8

**Caution:** Do not overinflate! Excessive pressure may rupture the system. Stop inflation when resistance is felt or the larynx is observed to close. Always release the pressure when finished training by opening the tubing port at the stopcock.

**Note:** Our Advanced Airway Management Trainer will fit on all models of **CPARLENE**<sup>®</sup> and non-modular manikins.

## Cleaning and Maintenance

To clean the Advanced Airway Management Trainer, you will first need to remove the head from the manikin. Next, take the trainer to an area with a sink and open counter space. Using the three red caps supplied, plug off the esophagus and the bronchi. **(See figure 8.)**

Stabilize the head on the counter face-up (towels work well for this), with the plugged tubes hanging over the sink. Carefully pour warm soapy water (a mild dish soap works best) into the mouth until the water level reaches halfway up the tongue. Now tilt the head back and bring the neck up 3" off the countertop.



Figure 9

Continue filling until the water level covers the tongue. At this point, take a small soft brush and gently scrub the inside of the mouth (a small toothbrush works well for this). Cotton swabs can be used to scrub inside the nostrils. When done, pull the plug from the esophagus and drain the water into the sink. Now pick the head up, hold it in a vertical position, and pull the plugs from the trachea to completely empty the system. **(See figure 9.)**

### Rinsing the Airway

To rinse the airway, follow the same procedure using clean, warm tap water. Repeat this process until all the soap has been flushed from the system.

### Disinfecting the Airway

To disinfect, repeat the standard cleaning procedure, but this time use a bleach solution instead of soapy water, as specified by the

Centers for Disease Control. Fill the system with the solution until it reaches the corners of the mouth. Remember to start filling with the head flat, and finish with the neck slightly elevated to ensure that the solution completely fills all airway passages. Once completely filled with the bleach solution, allow the head to sit for at least 10 minutes. Drain as described earlier and repeat the rinsing process to flush out all of the bleach solution. Set the head aside and allow it to dry completely.

### Cleaning the Face

Normal soil can be removed from the face of the trainer with mild, soapy water. Nasco Cleaner (LF09919U) will remove stubborn stains. Simply apply Nasco Cleaner to the soiled area and wipe clean with a soft cloth or paper towel.

**Note:** Do not use Nasco Cleaner around the mouth and nostrils, as the residue of the cleaner could be toxic.

### Caution

Never place the trainer on any kind of printed paper or plastic. These materials will transfer indelible stains. Ballpoint pens will also make indelible stains. Do not use cosmetics.

### Supplies/Replacement Parts for the Advanced “Airway Larry” Airway Management Trainer

**LF03605U** Replacement Lungs

**LF03606U** Replacement Stomach

**LF03644U** Pump Spray Lubricant

**LF09919U** Nasco Cleaner

## Other Available *Life/form*® Simulators

- LF03616U** Child *CRiSis*™ Manikin
- LF03617U** Deluxe Child *CRiSis*™ Manikin with Arrhythmia Tutor
- LF03620U** PALS Update Kit
- LF03699U** “Airway Larry” Airway Management Trainer
- LF03713U** **CPARLENE**®
- LF03953U** *CRiSis*™ Manikin, Complete
- LF03955U** Deluxe *CRiSis*™ Manikin
- LF03956U** Deluxe “Plus” *CRiSis*™ Manikin
- LF03965U** Adult *CRiSis*™ Auscultation Manikin
- LF03966U** Adult *CRiSis*™ Auscultation Manikin with ECG Simulator



*Nasco*

901 Janesville Avenue, P.O. Box 901  
Fort Atkinson, Wisconsin 53538-0901  
1.800.558.9595  
eNasco.com • E-mail: [lifeform@eNasco.com](mailto:lifeform@eNasco.com)