

# Endless love

## To Guard The New Life Safe



# YC-200A/200B

## Infant Bed



- Can be used in conjunction with jaundice treatment products, provide 24h bedside jaundice treatment and care for parents.



### CONFIGURATION TABLE ( "✓" Represents for standard function, "/" Represents for no function.)

	YC-200A	YC-200B
Infant basin tilt function	✓	/
Rotation arm swing function	✓	✓
Infant basin rotation function	/	✓
Retractable sidewall function	✓	✓
Crib height adjustment function	✓	✓



## Ningbo David Medical Device Co.,Ltd

Add: Bldg No. 5 Ningbo Smart Park, #98 Chuangyuan Rd, Hi-tech Zone, Ningbo, China

Tel: 0086-574-87800008 / 87800007

Fax: 0086-574-87803714

Web: [www.nbdavid.com/en/](http://www.nbdavid.com/en/)

E-mail: [sales@nbdavid.com](mailto:sales@nbdavid.com)

Specifications are changed without notice for performance improvement.

Printed date: Jul. 2022 Version No.: 1


## Rooming-in, creating a safe sleeping environment with love

With the family as the unit to reduce the risk of cross-infection, improve the exclusive breastfeeding rate, and establish the initial stage of parent-child interaction, promote the baby's intellectual growth and physical development.

### ■ 360° Rotating Infant Basin YC-200B


The 360° all-round rotation design, combined with the swing arm to stretch left and right, allows the mother to take care of the baby in a more comfortable position at a close distance, greatly enhance the parent-child relationship.

Mother can independently adjust the rotation angle to help lock in any comfortable position.

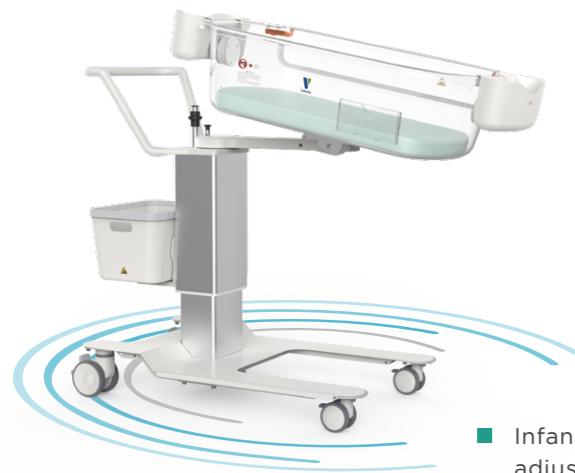

- Infant basin inclination angle is adjustable in three gears, the stable and safe structure design can avoid rollover even when tilted to the maximum angle.

YC-200A

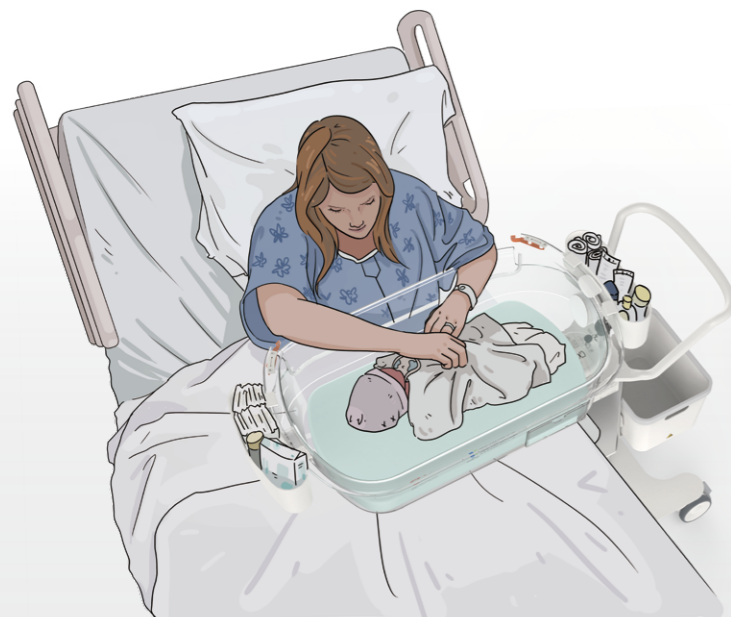


### ■ Recoverable side door panel

The use of scroll spring structure and small frictional resistance helps the mother to hold and care for the baby easily, which can minimize the noise generated, create a quiet and comfortable sleeping environment.

The design of door panel lock on both sides can prevent the baby from accidentally falling to the side and improve the safety level.

The door panel structure that can be easily disassembled, convenient for cleaning and disinfection without dead ends, which reduces the risk of infection.

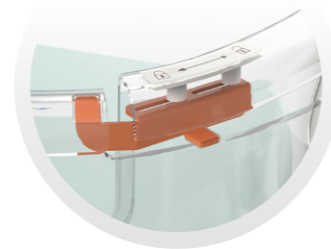

## Friendly design, provide intimate clinical care with love

Intimately meet the needs of mother to care for the baby during the postpartum recovery period. Even after caesarean operation, a mother can lie on the patient bed and follow the baby's movements at all times, leaving the love in zero distance.

- The transparent and visible baby basin made of PC material enables close interaction with the baby from a long distance even if the mother is lying on the bed.


- Pressure-relief anti-decubitus mattress is antibacterial, waterproof, soft and breathable, and protects the tender skin of newborns.



- The height of the bed is freely adjustable, and the ultra-high lifting stroke ( $27\text{cm} \pm 2\text{cm}$ ) can be adapted to any height of the hospital bed.

It can be easily locked at an ideal height to fully meet the needs of parent-child interaction and the nursing needs of medical staff of different heights.

- Bed frame in cold-rolled steel structure, which is sturdy and durable, protects the baby's nursing safety.

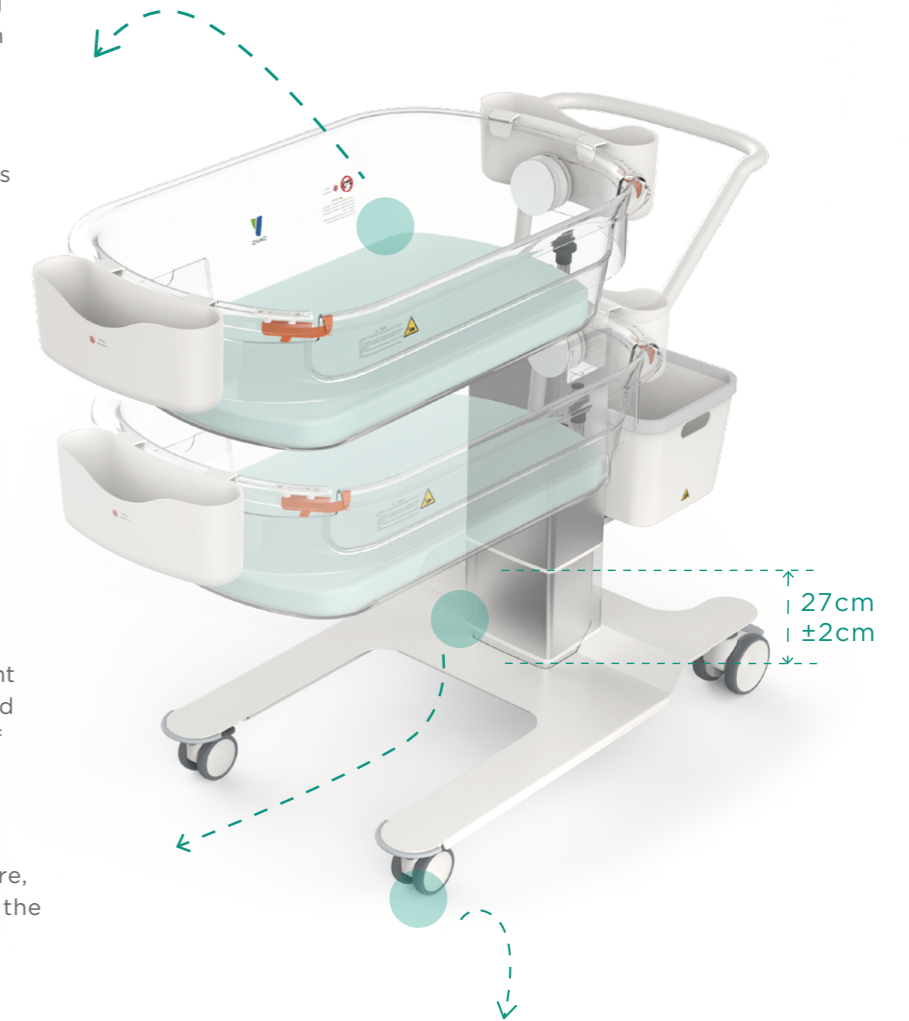

- Four universal locking casters are flexible and convenient to move, which can ensure the safety and stability during nursing.

- The large-capacity storage basket, together with the storage pockets on both sides of the bed, provides a larger storage space for maternal and child supplies and clinical care items.

Optional

A large-capacity storage basin is optional to meet more storage needs.