

Tuttnauer Dental

Bringing our
95 YEARS of Legacy
into your clinic



Quality Products ◆ Professional Excellence ◆ Reliable Services ◆ Time Saving Solutions

A complete instrument reprocessing solution for the patient-centric clinic.

Tailoring a smart solution for your clinic's steri-center with maximum automation and minimum time and cost.

A powerful combination of best practices and technology.





From **infection control** to **smart reprocessing**

The smart cleaning adaptable steri-center

Guiding you throughout the entire
process of instrument preparation.





Soaking & Washing

Washing the contaminated equipment

When packaging the equipment

Packaging & Sealing

Maintain sterility and prevent recontamination

Processing & Sterilization

Robust autoclave devices

Documentation & Traceability

Quality assurance for the contemporary clinic

Infection Control Indicators

Providing process quality throughout all the reprocessing steps

Tuttnauer Dental Premium Solution

Start

Soaking & Washing



Packaging & Sealing



Ultrasonic Cleaner



This quiet yet highly effective cleaner uses soundwaves to rid instruments of hard deposits.



Sealing Machine



To ensure effective, quality and long-lasting sterilization until the instruments are reused.






TIVA2 Washer-Disinfector



Simplify your work and ensure the safety of your patients and staff.

Automate your process, leave manual washing behind

- Fully automated, user-friendly operating and monitoring
- QA documentation, traceability and transparency
- Dozens of programs to choose from

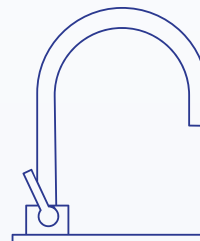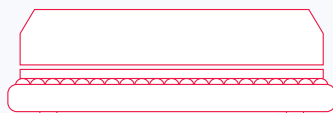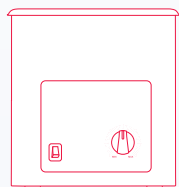

For our full range of models, please visit our website or open this QR code



Processing & Sterilization



T-Edge 10 Autoclave



Maximum efficiency. Minimum cost.

- Fully automated sterilization at the click of a button
- Chic design, excellent results, that fits your space limitation
- Simple-to-use interface
- Compliance with strict international standards
- Remote controlling and monitoring*

Comprehensive Set of Indicators

Details on the next page >>

Chamber Clean




Details on the next page >>

Weekly maintenance of the chamber.



Water Distiller

Small enough to fit on a counter; Smart enough to distill water for sterilization.

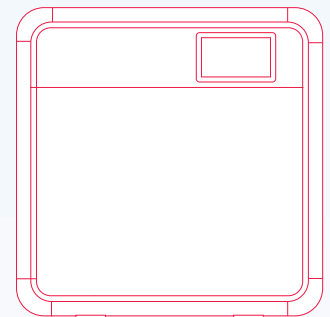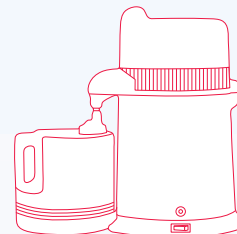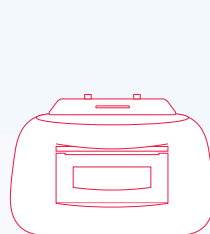

Documentation & Traceability

- Comprehensive set of QA tools
- Online traceability and documentation
- Biological and chemical indicators
- The only sterilization processor with transparent doors
- Receive barcodes with sterilization data (including expiration date).

*Remote monitoring capabilities will be available in the near future. Free upgrade will be offered under company terms.

Tuttnauer – Your Infection Control Partner

We deliver complete lines of indicators for ensuring comprehensive hygiene, cleaning, disinfecting, processing and sterilization.

Time saving: Highly effective conventional, rapid and super rapid indicators

Reliable performance: For the safety of your team and your patients

Easy to use: Simple to implement, read and record

Non-toxic materials: Protecting the environment and your staff

Extended shelf-life: Reducing waste and unnecessary expenses

Standards compliant: ISO 13485:2016 / NS-EN ISO 13485:2016 Quality Management System Standard, granted by the DNV International Certification Body



Biological Indicators



STEAM
WTL198-0058



STEAM
WTL198-0067



FLUORESCENCE
INCUBATOR
WTL198-0057



STEAM
WTL198-0056



INCUBATOR
PLZ198-0005

Fluorescence Super and Ultra Rapid readout

Conventional 24-hour readout

Chemical Indicators



WTL198-0053



WTL198-0055



WTL198-0068

Steam chemical
indicator tape

Steam chemical
indicator

Steam chemical
indicator

Washing & Disinfection Indicators



WTL198-0062



WTL198-0061



WTL198-0065



WTL198-0063



Pre-Sterilization Tests

Bowie-Dick and Helix

Type 2 Chemical Indicators

Tuttnauer offers a wide range of infection control products including kits to meet regulatory requirements and validate the dental reprocessing quality.



WTL198-0159

Helix-PCD + Bowie-Dick simulation test strips kit

- For monitoring vacuum-assisted steam sterilizers
- Autoclave release for hollow loads in large and small vacuum-assisted steam sterilizers
- Conditions: Bowie-Dick cycle, 3.5 minutes at 134 °C



WTL198-0161

Helix-PCD + chemical indicator test strips kit

- For safe release of hollow loads
- Process Challenge Device for Steam sterilization process control
- Conditions: 18 mins at 134 °C



WTL198-0158

Bowie-Dick test pack chemical indicator

- For detecting inadequate air removal and steam penetration
- Test pack for vacuum-assisted steam sterilizers
- Conditions: 16.9 mins at 121 °C; 4 mins at 132 °C; 3.5 minutes at 134 °C



WTL198-0160

Helix-PCD + chemical indicator test strips kit

- For safe release of hollow loads
- Process Challenge Device for steam sterilization processes control
- Conditions: 3.5 mins at 134 °C

Cost-effective and reliable autoclave monitoring
before and during sterilization cycles

Clean and Shine your Autoclave Chamber with our maintenance products

Chamber Brite

Debris free tablets to clean
and shine your Autoclave chamber.

For Class S & N Autoclaves



Chamber Clean

Hassle free tablets to clean
and shine your T-Edge chamber.

For T-Edge Class B Autoclaves



Your clinic's infection control success is backed by 95 years of experience of global leadership

Renowned global leaders in infection control solutions with more than

95 years experience

Products distributed to clinics, hospitals and research institutions in more than

140 countries

International Sales and Marketing
E-mail: info@tuttnauer.com
www.tuttnauer.com

Tuttnauer Europe b.v.
Hoeksteen 11, 4815 PR
PO Box 7191, 4800 GD Breda
The Netherlands
Tel: +31 765 423 510
Fax: +31 765 423 540
E-mail: info@tuttnauer.nl

Tuttnauer USA Co.
25 Power Drive,
Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850
Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com

Distributed by: